


Installation Instructions

"Series V" Control Panel

For King Air Compressors & Centrifugal Vacuums

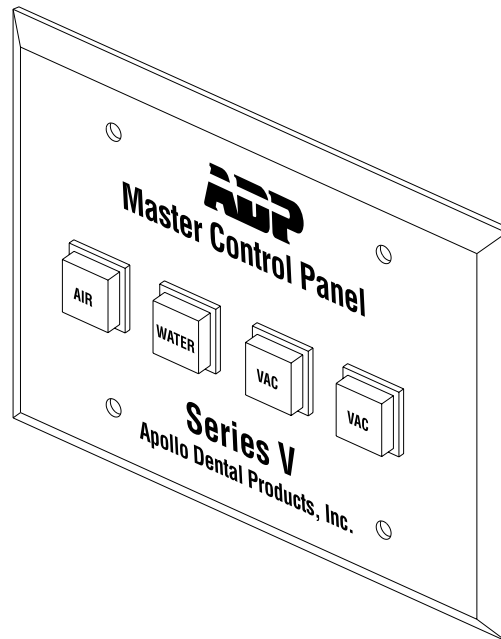

Description


The Apollo Dental Products' Control Panel are designed to provide a convenient means of controlling dental equipment from inside the office area. Features include lighted switches and an attractive brushed stainless steel mounting plate.

Master Control Panels include 50 feet of 18-3 bell wire with each switch, connecting fittings, and mounting fasteners.

Although the panel can be installed directly in drywall, many local codes require the use of a junction box. The junction box (**not provided**) must be a standard three gang box measuring 2 1/2" deep by 4" high by 5 3/4" wide. With either mounting method, the panel requires a standard three wire **24 VAC** control circuit from each piece of equipment being controlled. ADP also manufactures a complete line of remote contactor boxes for **24 VAC** control.

Instructions For Junction Box Mounting

1. Prior to installation of sheet rock, mount the three gang junction box to a stud as illustrated at right. The box should be placed so that it's front surface is set back 1/2" to 5/8".
2. Run one 18-3 bell wire from each piece of equipment up to the installed junction box.
3. Cut and install sheet rock.
4. Connectors are provided to complete the wiring of the control panel by connecting wires to corresponding colored wires from the Apollo low voltage contactor box(es).
5. Position the panel, then insert and tighten screws.

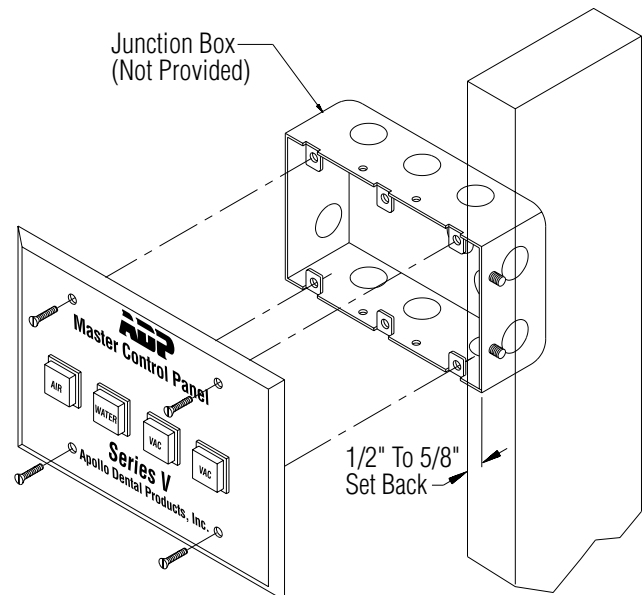


Box Mounted Control Panel Diagram

Instructions For Drywall Mounting


Warning: All Control Panels must be installed according to local codes.

1. Place the provided template on wall and mark the position of the four 3/16" diameter holes and the switch area.
2. Cut and remove the area from the drywall.
CAUTION: Do not cut hole larger than shown on the template.
3. Drill 3/16" holes.
4. Insert the provided anchors into the 3/16" holes.
5. Run one 18-3 bell wire from each piece of equipment to the area the panel is to be mounted.
6. Use the provided connectors to complete the wiring of the control panel. Attach each colored wire to the corresponding colored wire from the ADP contactor box(es).
7. Position the panel and insert screws. After the screw has contacted the mounting plate, use the chart below to find the additional number of turns needed for proper expansion of the anchors.

Drywall Thickness	Additional Turns After Contact
1/2"	5 - 7 Turns
3/4"	3 - 4 Turns
1"	1 - 2 Turns

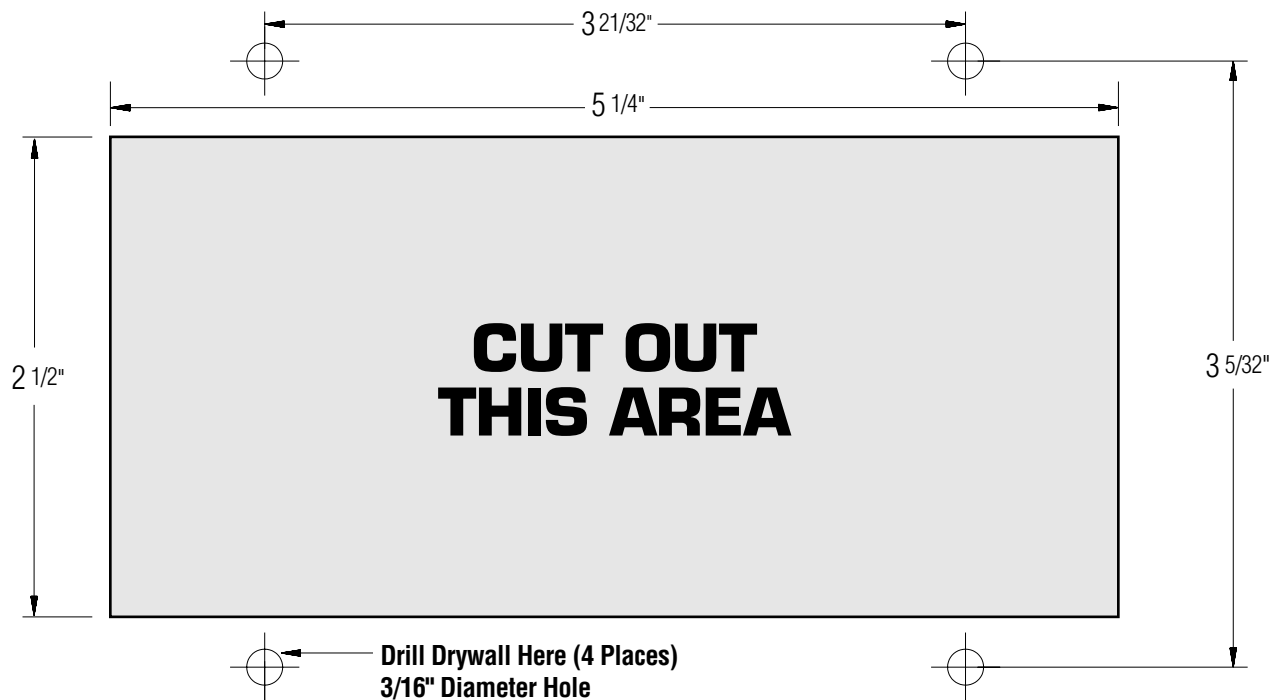

Drywall Mounting Template

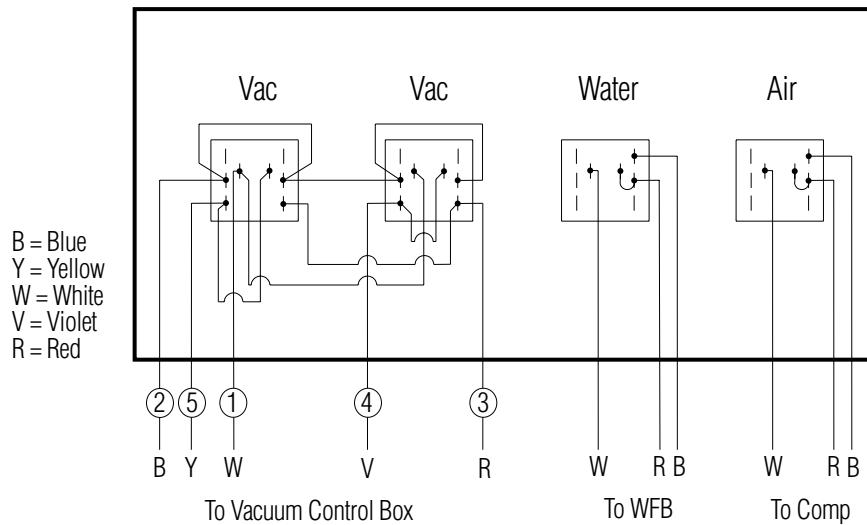

Note: Viewed from the back side

Twin Schematic (4 Switch)

Trouble-Shooting Guide

Problem: *Dental equipment does not start when switch is engaged.*

Cause A: Broken or loose low voltage wire.

Corrective Action: 1. Check all wire connections.

Cause B: Defective switch.

Corrective Action: 1. Connect red and blue wires together at the conntactor box, thus bypassing the switch.

Problem: Dental equipment starts when switch is engaged but switch does not illuminate.

Cause A: Broken or loose wire.

Corrective Action: 1. Check all wire connections.

Cause B: Defective bulb.

Corrective Action: 1. Relamp switch by removing switch lens, then pulling bulb ejector lever located on inner lower edge of switch. Press in new bulb (#SEA95225).

Warranty Information: 90 Days

All ADP accessories are thoroughly inspected and tested in accordance with rigid specifications and standards. Our accessories are guaranteed against any defective material and workmanship from the date of shipment; provided that the installation, operation, and maintenance is done in accordance with ADP procedures as outlined in our Installation and Maintenance Guides. Warranty cards must be returned to ADP within ten days of installation to effect warranty. No other warranties or guarantees, expressed or implied are made.

ADP's obligation under the warranty is to provide parts for the repair or, at its option, to provide the replacement product (excluding labor). All special, incidental and/or consequential damages are excluded. We will not issue credit for accessories without first attempting to correct the problem in the field. Written notice of breach of warranty must be given to ADP within the warranty period. The warranty does not cover damage resulting from improper installation or maintenance, accident or misuse. The warranty does not cover damage resulting from the use of cleaning, disinfecting or sterilizing chemicals and processes. The warranty does not cover vacuum failures due to hard water deposits. Failure to follow instructions provided in ADP's Installation and Maintenance Guides may void the warranty.